Ambu[®] Perfit ACE[®] & Mini Perfit ACE[®]



Key Benefits

- Stores completely flat
- CT, MR and X-ray compatible*
- Standard size markings on both Adult and Mini
- Chin flip feature to ease airway management
- Tube hooks are located on the collar body to secure the oxygen tubes
- Facilitates intubation
- Disposable
- When shipped, the default position of the the Adult Ace collar is in the Neckless position, and the default position of the Mini Ace is the Small Pediatric position

Ambu[®] Perfit ACE[®] & Mini Perfit ACE[®]

28 sizes in two collars

The Ambu Perfit ACE is a unique adjustable collar. Combined with the Mini Perfit ACE these two collars offer an incredible 28 sizes ensuring that patients can now have the most effective and accurate fit ever.

The collar is packaged completely flat and has a unique chin-flip feature to ease airway management. Ambu Perfit ACE is a disposable adult adjustable collar for immobilizing the cervical column.

Ambu Perfit Ace is an adjustable collar designed to fit any adult with 16 different settings ranging from Neckless to Tall sizes. The Ambu Perfit ACE incorporates ratchet latches that allow the rescuer to size the collar exactly to the patient's neck size. Two locking buttons hold the collar in place when adjusted.

Ambu Mini Perfit ACE is an adjustable collar designed to fit small pediatric, pediatric and small adult, with 12 different settings to provide an accurate fit.

* Distortion or shadows may appear on the image depending on the settings of the imaging equipment



Ambu Perfit ACE

Adjustable sizes	Neckless, Short, Regular and Tall
Number of settings	16 adjustments
Dimensions	554 x 215 x 14 mm (21.8" x 8.5" x 0.6")
Weight	200 g (7 ounces)
Color	White/light grey
Ambu Mini Perfit ACE	
Adjustable sizes	Small pediatric, Pediatric and Small Adult
Number of settings	12 adjustments
Dimensions	447 x 188 x 14 mm (17.6" x 7.4" x 0.6")
Weight	100 g (3.5 ounces)
Color	Yellow/light grey



Ambu[®] Perfit ACE[®]


Material

Product	Material
Ambu Perfit ACE	
White/light grey plastic	Polypropylene
Foam	Polyethylene
Ambu Mini Perfit Ace	
Yellow plastic	Polypropylene
Foam	Polyethylene

Ordering specifications

Item no.	Description
002 810 001	Perfit ACE Extrication Collar, 20/Case
000 281 107	Mini Perfit ACE Extrication Collar, 30/Case

Ambu® Mini Perfit ACE®

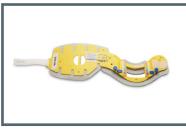

Ambu[®] Mini Perfit ACE[®]

Ambu A/S Baltorpbakken 13 DK - 2750 Ballerup Denmark Tel. +45 7225 2000 ambu.com