

Ambu® AuraStraight™



Key Benefits

- Reinforced tip resists folding over during insertion and plugs the upper esophageal sphincter
- Cuff and airway tube molded as single unit for extra safety - no separation
- The surface has an EasyGlide texture to ease insertion
- Extra soft cuff ensures the best possible seal with least possible mucosal pressure
- Pilot balloon identifies mask size and provides precise tactile indication of degree of inflation
- Convenient depth marks for monitoring correct position of the mask
- Packaged sterile and ready for use
- Color coded pouch with directions for use
- Available in 8 sizes fitting a wide range of patients
- MR safe
- Does not contain natural rubber latex

Ambu® AuraStraight™









Ambu AuraStraight is a disposable laryngeal mask designed for use both in anesthesia and emergency medicine. It provides users with an even greater choice and is an alternative to the already popular AuraOnce curved laryngeal mask.


Unlike other laryngeal masks, AuraStraight features a unique one-piece-mold design, which ensures not only consistency, but also that there are no rough edges on the airway tube or cuff.

Furthermore, the AuraStraight provides an extra soft, thin cuff. Thus, the seal more readily conforms to the shape of the airway with the minimal amount of pressure. Additionally, the reinforced tip and EasyGlide surface aids the ease of insertion and extraction.

Ambu AuraStraight comes in a full range of 8 sizes from neonate to large adult.

Specifications

	Infant		Pediatric		Adult			
Mask size	#1	#1 1/2	#2	#2½	#3	#4	#5	#6
Packaging colour code								
Patent weight	<5 kg	5-10 kg	10-20 kg	20-30 kg	30-50 kg	50-70 kg	70-100 kg	>100 kg
Maximum cuff inflation volume	4 ml	7 ml	10 ml	14 ml	20 ml	30 ml	40 ml	50 ml
Maximum intra-cuff pressure	60 cm H ₂ O							
Airway connector	15 mm ISO 5356-1							
Min. I.D. Tube	5.7 mm	6.5 mm	7.4 mm	8.9 mm	9.8 mm	10.7 mm	11.6 mm	12.1 mm
Max. O.D. Tube	9.6 mm	11 mm	12.2 mm	14.4 mm	16.1 mm	17.6 mm	19.3 mm	19.5 mm
Inflation valve	Luer cone ISO 594-1							
Appropriate storage temperature	+10° C - +25° C							
Weight	9.5 g	13.5 g	17.5 g	26.6 g	37.1 g	46 g	63.4 g	69.7 g
Internal volume of ventilatory pathway	3.5 ml	5.4 ml	6.7 ml	9.7 ml	13 ml	15 ml	20 ml	22 ml
Pressure drop	<1.5 cm H ₂ O at 15 l/min	<0.8 cm H ₂ O at 15 l/min	<1.9 cm H ₂ O at 30 l/min	<1.0 cm H ₂ O at 30 l/min	<2.2 cm H ₂ O at 60 l/min	<1.8 cm H ₂ O at 60 l/min	<1.2 cm H ₂ O at 60 l/min	<1.2 cm H ₂ O at 60 l/min
Min. interdental gap	12.8 mm	15.1 mm	17.1 mm	20 mm	23 mm	25.4 mm	27.4 mm	28.7 mm
Internal Pathway	10.9 cm	11.6 cm	13 cm	14.4 cm	16 mm	17.1 mm	18.8 mm	20.4 cm
Specialized use								
Max. FOB size	2.7 mm	2.7 mm	3.9 mm	3.9 mm	5.4 mm	5.4 mm	5.4 mm	5.4 mm
Max. ETT size	3.0 mm	3.5 mm	4.0 mm	5.0 mm	5.5 mm	6.0 mm	6.5 mm	6.5 mm

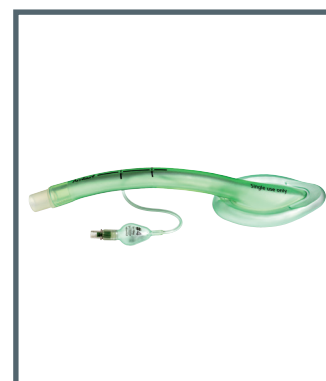

AuraStraight



AuraStraight with reinforced tip

Materials

Product		Material
Airway tube		PVC
Cuff		PVC
Airway connector		Polypropylen (PP)
Pilot Balloon		PVC
Inflation valve		PVC/Silicone
Pilot tube		PVC
Packaging	Quantity	Material
Pouch	1	Tyvek®/Polyethylene Terephthalate (PET)/Polyethylene (PE)
Box	10	Cardboard

Ordering specifications

Item no.	Size	Type
324 100 000	1	Disposable Laryngeal Mask - Straight
324 150 000	1 1/2	Disposable Laryngeal Mask - Straight
324 200 000	2	Disposable Laryngeal Mask - Straight
324 250 000	2 1/2	Disposable Laryngeal Mask - Straight
324 300 000	3	Disposable Laryngeal Mask - Straight
324 400 000	4	Disposable Laryngeal Mask - Straight
324 500 000	5	Disposable Laryngeal Mask - Straight
324 600 000	6	Disposable Laryngeal Mask - Straight



AuraStraight in colour-coded pouch with instructions for use

Ambu Australia Pty Ltd

Customer Services

Tel.: 1300 233 118
Fax: 1300 156 526
aussales@ambu.com

Ambu Asia (Sales office)

Customer Services

Tel.: +604 252 9100
Fax: +604 644 8932
csasia@ambu.com